

**KAMIYA BIOMEDICAL COMPANY**

# Dog Microalbumin ELISA

**For the quantitative determination of microalbumin  
in dog biological fluids**

**Cat. No. KT-389**

**For Research Use Only.**

**PRODUCT INFORMATION****Dog Microalbumin ELISA**  
**Cat. No. KT-389****INTENDED USE**

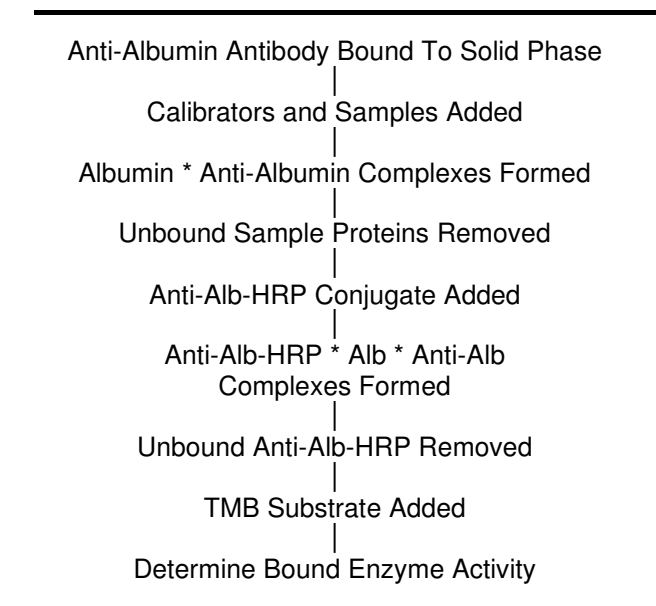
The Dog Microalbumin ELISA is a highly sensitive two-site enzyme-linked immunoassay (ELISA) for the quantitative determination of microalbumin in dog biological fluids. For research use only.

**INTRODUCTION**

Albumin (Alb) is an amazing polyfunctional protein contributing to homeostasis through mechanisms of hemodynamics, transport and nutrition. Albumin is found both intra- and extra-vascularly in all mammals and many lower vertebrates. It is a molecule of about 67,000 daltons, synthesized by the liver. Normally only very trace amounts of albumin escape reabsorption by kidney glomeruli and are excreted into the urine. Many occult diseases can cause kidney damage that may result in excessive amounts of serum proteins, including albumin, to be excreted by the kidney and into the urine. This ELISA can be used to measure albumin in serum, tissue extracts and other biological fluids.

**PRINCIPLE**

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the albumin present in the sample reacts with the anti-Alb antibody, which has been adsorbed to the surface of polystyrene microtiter wells. After the removal of unbound proteins by washing, anti-Alb antibody conjugated with horseradish peroxidase (HRP) is added. This HRP-conjugated antibody forms a complex with the previously bound albumin. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme is proportional to the concentration of albumin in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of albumin in the test sample. The quantity of albumin in the test sample can be interpolated from the calibration curve constructed from the calibrators and corrected for sample dilution.

**Figure 1.****COMPONENTS**

1. Diluent Concentrate  
One bottle containing 50 mL of a 5X concentrated diluent running buffer.
2. Wash Solution Concentrate  
One bottle containing 50 mL of a 20X concentrated wash solution.

3. Enzyme-Antibody Conjugate Concentrate  
One vial containing 150  $\mu$ L of a 100X concentrated affinity-purified anti-dog albumin antibody conjugated with HRP in a stabilizing buffer.
4. TMB Substrate Solution  
One vial containing 12 mL of TMB and hydrogen peroxide in citric acid buffer at pH 3.3.
5. Stop Solution  
One vial containing 12 mL of 0.3 M sulfuric acid.  
**WARNING: Avoid contact with skin.**
6. Microtiter Plate  
Twelve removable eight-well strips in well holder frame. Wells are coated with affinity-purified anti-dog albumin.
7. Dog Microalbumin Calibrator  
One vial containing a lyophilized Dog Microalbumin Calibrator.
8. Positive Control  
One vial containing 50  $\mu$ L of serum with 0.1% sodium azide. See the control certificate for the concentration.

### MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes (2  $\mu$ L to 200  $\mu$ L) for making and dispensing dilutions
- Test tubes
- Microplate washer/aspirator
- Distilled or de-ionized H<sub>2</sub>O
- Microplate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer

### PRECAUTIONS

1. Read the instructions carefully before beginning the assay.
2. This kit is for research use only.
3. Great care has been taken to ensure the quality and reliability of this product. However, it is possible that in certain cases, unusual results may be obtained due to high levels of interfering factors.
4. Preservatives
  - Positive Control contains 0.1% sodium azide.
5. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.
6. Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.
7. Other precautions:
  - Do not interchange kit components from different lots.
  - Do not use kit components beyond the expiration date.
  - Protect reagents from direct sunlight.
  - Do not pipette by mouth.
  - Do not eat, drink, smoke or apply cosmetics where reagents are used.
  - Avoid all contact with the reagents by using gloves.
  - Stop solution contains diluted sulfuric acid. Irritation to eyes and skin is possible. Flush with water after contact.

### REAGENT PREPARATION

1. Diluent Concentrate  
The Diluent solution supplied is a 5X concentrate and must be diluted 1:5 with distilled or de-ionized water.
2. Wash Solution Concentrate  
The Wash Solution supplied is a 20X concentrate and must be diluted 1:20 with distilled or de-ionized water. Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.
3. Enzyme-Antibody Conjugate Concentrate

Calculate the required amount of working conjugate solution for each microtiter plate test strip by adding 10  $\mu$ L Enzyme-Antibody Conjugate to 990  $\mu$ L of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.

4. TMB Substrate Solution  
Ready to use as supplied.
5. Stop Solution  
Ready to use as supplied.
6. Microtiter Plate  
Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.
7. Dog Microalbumin Calibrator  
Add 1.0 mL of distilled or de-ionized water to the lyophilized Dog Microalbumin Calibrator and mix gently until dissolved. The calibrator is now at a concentration of 233.333  $\mu$ g/mL (**the reconstituted calibrator should be aliquoted and frozen if future use is intended**). Dog Microalbumin Calibrators need to be prepared immediately prior to use (see chart below). Mix well between each step. Avoid foaming.

Calibrator	Concentration (ng/mL)	Calibrator Volume added to 1X Diluent	Volume of 1X Diluent
A	2,333.33	5 $\mu$ L Dog Microalbumin Calibrator	495 $\mu$ L
6	400	150 $\mu$ L Calibrator A	725 $\mu$ L
5	200	0.3 mL Calibrator 6	0.3 mL
4	100	0.3 mL Calibrator 5	0.3 mL
3	50	0.3 mL Calibrator 4	0.3 mL
2	25	0.3 mL Calibrator 3	0.3 mL
1	12.5	0.3 mL Calibrator 2	0.3 mL
0	0		0.6 mL

8. Positive Control  
The concentration and recommended dilution are provided on the control certificate. Before use, briefly centrifuge the Positive Control to allow all of the liquid to collect in the bottom of the vial.

## STORAGE AND STABILITY

1. Complete Kit  
The expiration date for the kit is stated on the outer label. The recommended storage temperature is 4°C. **Note: See long term storage recommendations below for the Dog Microalbumin Calibrator and Positive Control.**
2. Diluent  
The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
3. Wash Solution  
The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
4. Enzyme-Antibody Conjugate  
Undiluted anti-Alb-HRP conjugate should be stored at 4°C and **diluted immediately prior to use**. The working conjugate solution is stable for up to 1 hour if stored in the dark.
5. TMB Substrate Solution  
The TMB Substrate Solution should be stored at 4°C and is stable until the expiration date.
6. Stop Solution  
The Stop Solution should be stored at 4°C and is stable until the expiration date.

7. **Microtiter Plate**  
Anti-dog albumin coated wells are stable until the expiration date, and should be stored at 4°C in the sealed foil pouch with desiccant pack.
8. **Dog Microalbumin Calibrator**  
**The lyophilized Dog Microalbumin Calibrator should be stored at 4°C or frozen until reconstituted. The reconstituted calibrator should be aliquoted and stored frozen. Avoid multiple freeze/thaw cycles.** The working calibrator solutions should be prepared immediately prior to use and are stable for up to 8 hours.
9. **Positive Control**  
For storage longer than 7 days keep frozen until the expiration date. Storage less than 7 days can be at 4°C. Avoid multiple freeze/thaw cycles.

## INDICATIONS OF INSTABILITY

If the test is performing correctly, the results observed with the calibrator solutions should be within 20% of the expected values.

## SPECIMEN COLLECTION AND HANDLING

Blood should be collected by venipuncture and the serum separated from the cells, after clot formation, by centrifugation. For plasma samples, blood should be collected into a container with an anticoagulant and then centrifuged. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. For urine samples, collect using standard techniques. The samples should be centrifuged to remove any possible debris. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freezing/thawing.

For any sample that might contain pathogens, care must be taken to prevent contact with open wounds. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.

## ASSAY PROTOCOL

### Dilution of Samples

The assay for quantification of albumin in urine, serum or plasma requires that each sample be diluted before use. A 1:1,000 dilution is appropriate for most urine samples while serum or plasma samples may need to be diluted 1:200,000. For absolute quantification of samples that yield results outside the range of the calibration curve, a lesser or greater dilution might be required. If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.

To prepare a 1:1,000 dilution of sample, transfer 5 µL of sample to 495 µL of 1X Diluent. This gives you a 1:100 dilution. Then, transfer 40 µL of your 1:100 dilution to 360 µL of 1X Diluent. You now have a 1:1,000 dilution. To prepare a 1:200,000 dilution, transfer 2 µL of your 1:1,000 dilution to 398 µL of 1X Diluent. You now have a 1:200,000 dilution. Make sure to mix thoroughly at each stage.

### Procedure

1. **Bring all reagents to RT before use.**
2. Pipette 100 µL of
  - Calibrator 0 (0.0 ng/mL) in duplicate
  - Calibrator 1 (12.5 ng/mL) in duplicate
  - Calibrator 2 (25 ng/mL) in duplicate
  - Calibrator 3 (50 ng/mL) in duplicate
  - Calibrator 4 (100 ng/mL) in duplicate
  - Calibrator 5 (200 ng/mL) in duplicate
  - Calibrator 6 (400 ng/mL) in duplicate
3. Pipette 100 µL of Positive Control (in duplicate) into pre-designated wells.
4. Pipette 100 µL of sample (in duplicate) into pre-designated wells.

5. Incubate the Microtiter Plate at 22°C (RT) for sixty (60 ± 2) minutes. Keep plate covered and level during incubation.
6. Following incubation, aspirate the contents of the wells.
7. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with diluted Wash Solution, invert the plate and pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual Wash Solution. Repeat three times for a total of four washes.
8. Pipette 100 µL of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at 22°C (RT) for thirty (30 ± 2) minutes. Keep plate covered in the dark and level during incubation.
9. Wash and blot the wells as described in Steps 6 and 7.
10. Pipette 100 µL of TMB Substrate Solution into each well.
11. Incubate in the dark at RT for precisely ten (10) minutes.
12. After ten minutes, add 100 µL of Stop Solution to each well.
13. Determine the absorbance at 450 nm of the contents of each well. Calibrate the plate reader to manufacturer's specifications.

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

## RESULTS

1. Subtract the average background value from the test values for each sample.
2. Using the results observed for the calibrators construct a calibration curve. The appropriate curve fit is that of a four-parameter logistics curve, although a second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from calibration curve. Correct for sample dilution factor to arrive at microalbumin concentration in original sample.

## PERFORMANCE CHARACTERISTICS

In accord with good laboratory practice, the assays for specific microalbumin require meticulous quality control. Each laboratory should use routine quality control procedures to establish inter- and intra-assay precision and performance characteristics.

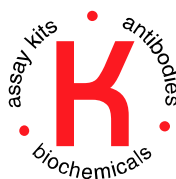
## LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or de-ionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.

## **FOR RESEARCH USE ONLY**

### **KAMIYA BIOMEDICAL COMPANY**

12779 Gateway Drive, Seattle, WA 98168  
Tel: (206) 575-8068 Fax: (206) 575-8094  
Email: LifeScience@k-assay.com  
www.k-assay.com



**KAMIYA BIOMEDICAL COMPANY**

# Dog Microalbumin ELISA

**For the quantitative determination of microalbumin  
in dog biological fluids**

**Cat. No. KT-389**

**For Research Use Only.**

**PRODUCT INFORMATION****Dog Microalbumin ELISA**  
**Cat. No. KT-389****INTENDED USE**

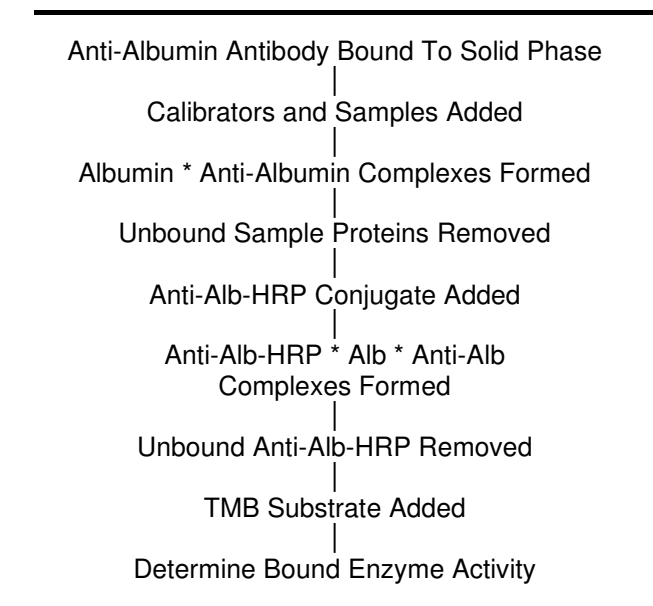
The Dog Microalbumin ELISA is a highly sensitive two-site enzyme-linked immunoassay (ELISA) for the quantitative determination of microalbumin in dog biological fluids. For research use only.

**INTRODUCTION**

Albumin (Alb) is an amazing polyfunctional protein contributing to homeostasis through mechanisms of hemodynamics, transport and nutrition. Albumin is found both intra- and extra-vascularly in all mammals and many lower vertebrates. It is a molecule of about 67,000 daltons, synthesized by the liver. Normally only very trace amounts of albumin escape reabsorption by kidney glomeruli and are excreted into the urine. Many occult diseases can cause kidney damage that may result in excessive amounts of serum proteins, including albumin, to be excreted by the kidney and into the urine. This ELISA can be used to measure albumin in serum, tissue extracts and other biological fluids.

**PRINCIPLE**

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the albumin present in the sample reacts with the anti-Alb antibody, which has been adsorbed to the surface of polystyrene microtiter wells. After the removal of unbound proteins by washing, anti-Alb antibody conjugated with horseradish peroxidase (HRP) is added. This HRP-conjugated antibody forms a complex with the previously bound albumin. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme is proportional to the concentration of albumin in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of albumin in the test sample. The quantity of albumin in the test sample can be interpolated from the calibration curve constructed from the calibrators and corrected for sample dilution.

**Figure 1.****COMPONENTS**

1. Diluent Concentrate  
One bottle containing 50 mL of a 5X concentrated diluent running buffer.
2. Wash Solution Concentrate  
One bottle containing 50 mL of a 20X concentrated wash solution.



3. Enzyme-Antibody Conjugate Concentrate  
One vial containing 150  $\mu$ L of a 100X concentrated affinity-purified anti-dog albumin antibody conjugated with HRP in a stabilizing buffer.
4. TMB Substrate Solution  
One vial containing 12 mL of TMB and hydrogen peroxide in citric acid buffer at pH 3.3.
5. Stop Solution  
One vial containing 12 mL of 0.3 M sulfuric acid.  
**WARNING: Avoid contact with skin.**
6. Microtiter Plate  
Twelve removable eight-well strips in well holder frame. Wells are coated with affinity-purified anti-dog albumin.
7. Dog Microalbumin Calibrator  
One vial containing a lyophilized Dog Microalbumin Calibrator.
8. Positive Control  
One vial containing 50  $\mu$ L of serum with 0.1% sodium azide. See the control certificate for the concentration.

### MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes (2  $\mu$ L to 200  $\mu$ L) for making and dispensing dilutions
- Test tubes
- Microplate washer/aspirator
- Distilled or de-ionized H<sub>2</sub>O
- Microplate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer

### PRECAUTIONS

1. Read the instructions carefully before beginning the assay.
2. This kit is for research use only.
3. Great care has been taken to ensure the quality and reliability of this product. However, it is possible that in certain cases, unusual results may be obtained due to high levels of interfering factors.
4. Preservatives
  - Positive Control contains 0.1% sodium azide.
5. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.
6. Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.
7. Other precautions:
  - Do not interchange kit components from different lots.
  - Do not use kit components beyond the expiration date.
  - Protect reagents from direct sunlight.
  - Do not pipette by mouth.
  - Do not eat, drink, smoke or apply cosmetics where reagents are used.
  - Avoid all contact with the reagents by using gloves.
  - Stop solution contains diluted sulfuric acid. Irritation to eyes and skin is possible. Flush with water after contact.

### REAGENT PREPARATION

1. Diluent Concentrate  
The Diluent solution supplied is a 5X concentrate and must be diluted 1:5 with distilled or de-ionized water.
2. Wash Solution Concentrate  
The Wash Solution supplied is a 20X concentrate and must be diluted 1:20 with distilled or de-ionized water. Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.
3. Enzyme-Antibody Conjugate Concentrate

Calculate the required amount of working conjugate solution for each microtiter plate test strip by adding 10  $\mu$ L Enzyme-Antibody Conjugate to 990  $\mu$ L of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.

4. TMB Substrate Solution  
Ready to use as supplied.
5. Stop Solution  
Ready to use as supplied.
6. Microtiter Plate  
Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.
7. Dog Microalbumin Calibrator  
Add 1.0 mL of distilled or de-ionized water to the lyophilized Dog Microalbumin Calibrator and mix gently until dissolved. The calibrator is now at a concentration of 233.333  $\mu$ g/mL (**the reconstituted calibrator should be aliquoted and frozen if future use is intended**). Dog Microalbumin Calibrators need to be prepared immediately prior to use (see chart below). Mix well between each step. Avoid foaming.

Calibrator	Concentration (ng/mL)	Calibrator Volume added to 1X Diluent	Volume of 1X Diluent
A	2,333.33	5 $\mu$ L Dog Microalbumin Calibrator	495 $\mu$ L
6	400	150 $\mu$ L Calibrator A	725 $\mu$ L
5	200	0.3 mL Calibrator 6	0.3 mL
4	100	0.3 mL Calibrator 5	0.3 mL
3	50	0.3 mL Calibrator 4	0.3 mL
2	25	0.3 mL Calibrator 3	0.3 mL
1	12.5	0.3 mL Calibrator 2	0.3 mL
0	0		0.6 mL

8. Positive Control  
The concentration and recommended dilution are provided on the control certificate. Before use, briefly centrifuge the Positive Control to allow all of the liquid to collect in the bottom of the vial.

## STORAGE AND STABILITY

1. Complete Kit  
The expiration date for the kit is stated on the outer label. The recommended storage temperature is 4°C. **Note: See long term storage recommendations below for the Dog Microalbumin Calibrator and Positive Control.**
2. Diluent  
The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
3. Wash Solution  
The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.
4. Enzyme-Antibody Conjugate  
Undiluted anti-Alb-HRP conjugate should be stored at 4°C and **diluted immediately prior to use**. The working conjugate solution is stable for up to 1 hour if stored in the dark.
5. TMB Substrate Solution  
The TMB Substrate Solution should be stored at 4°C and is stable until the expiration date.
6. Stop Solution  
The Stop Solution should be stored at 4°C and is stable until the expiration date.

7. **Microtiter Plate**  
Anti-dog albumin coated wells are stable until the expiration date, and should be stored at 4°C in the sealed foil pouch with desiccant pack.
8. **Dog Microalbumin Calibrator**  
**The lyophilized Dog Microalbumin Calibrator should be stored at 4°C or frozen until reconstituted. The reconstituted calibrator should be aliquoted and stored frozen. Avoid multiple freeze/thaw cycles.** The working calibrator solutions should be prepared immediately prior to use and are stable for up to 8 hours.
9. **Positive Control**  
For storage longer than 7 days keep frozen until the expiration date. Storage less than 7 days can be at 4°C. Avoid multiple freeze/thaw cycles.

## INDICATIONS OF INSTABILITY

If the test is performing correctly, the results observed with the calibrator solutions should be within 20% of the expected values.

## SPECIMEN COLLECTION AND HANDLING

Blood should be collected by venipuncture and the serum separated from the cells, after clot formation, by centrifugation. For plasma samples, blood should be collected into a container with an anticoagulant and then centrifuged. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. For urine samples, collect using standard techniques. The samples should be centrifuged to remove any possible debris. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freezing/thawing.

For any sample that might contain pathogens, care must be taken to prevent contact with open wounds. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.

## ASSAY PROTOCOL

### Dilution of Samples

The assay for quantification of albumin in urine, serum or plasma requires that each sample be diluted before use. A 1:1,000 dilution is appropriate for most urine samples while serum or plasma samples may need to be diluted 1:200,000. For absolute quantification of samples that yield results outside the range of the calibration curve, a lesser or greater dilution might be required. If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.

To prepare a 1:1,000 dilution of sample, transfer 5 µL of sample to 495 µL of 1X Diluent. This gives you a 1:100 dilution. Then, transfer 40 µL of your 1:100 dilution to 360 µL of 1X Diluent. You now have a 1:1,000 dilution. To prepare a 1:200,000 dilution, transfer 2 µL of your 1:1,000 dilution to 398 µL of 1X Diluent. You now have a 1:200,000 dilution. Make sure to mix thoroughly at each stage.

### Procedure

1. **Bring all reagents to RT before use.**
2. Pipette 100 µL of
  - Calibrator 0 (0.0 ng/mL) in duplicate
  - Calibrator 1 (12.5 ng/mL) in duplicate
  - Calibrator 2 (25 ng/mL) in duplicate
  - Calibrator 3 (50 ng/mL) in duplicate
  - Calibrator 4 (100 ng/mL) in duplicate
  - Calibrator 5 (200 ng/mL) in duplicate
  - Calibrator 6 (400 ng/mL) in duplicate
3. Pipette 100 µL of Positive Control (in duplicate) into pre-designated wells.
4. Pipette 100 µL of sample (in duplicate) into pre-designated wells.

5. Incubate the Microtiter Plate at 22°C (RT) for sixty (60 ± 2) minutes. Keep plate covered and level during incubation.
6. Following incubation, aspirate the contents of the wells.
7. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with diluted Wash Solution, invert the plate and pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual Wash Solution. Repeat three times for a total of four washes.
8. Pipette 100 µL of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at 22°C (RT) for thirty (30 ± 2) minutes. Keep plate covered in the dark and level during incubation.
9. Wash and blot the wells as described in Steps 6 and 7.
10. Pipette 100 µL of TMB Substrate Solution into each well.
11. Incubate in the dark at RT for precisely ten (10) minutes.
12. After ten minutes, add 100 µL of Stop Solution to each well.
13. Determine the absorbance at 450 nm of the contents of each well. Calibrate the plate reader to manufacturer's specifications.

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

## RESULTS

1. Subtract the average background value from the test values for each sample.
2. Using the results observed for the calibrators construct a calibration curve. The appropriate curve fit is that of a four-parameter logistics curve, although a second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from calibration curve. Correct for sample dilution factor to arrive at microalbumin concentration in original sample.

## PERFORMANCE CHARACTERISTICS

In accord with good laboratory practice, the assays for specific microalbumin require meticulous quality control. Each laboratory should use routine quality control procedures to establish inter- and intra-assay precision and performance characteristics.

## LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or de-ionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.

## **FOR RESEARCH USE ONLY**

### **KAMIYA BIOMEDICAL COMPANY**

12779 Gateway Drive, Seattle, WA 98168  
Tel: (206) 575-8068 Fax: (206) 575-8094  
Email: LifeScience@k-assay.com  
www.k-assay.com