

KAMIYA BIOMEDICAL COMPANY

Rat Kidney Injury Molecule (KIM1) ELISA

**For the quantitative determination of KIM1
in Rat biological fluids**

Cat. No. KT-622

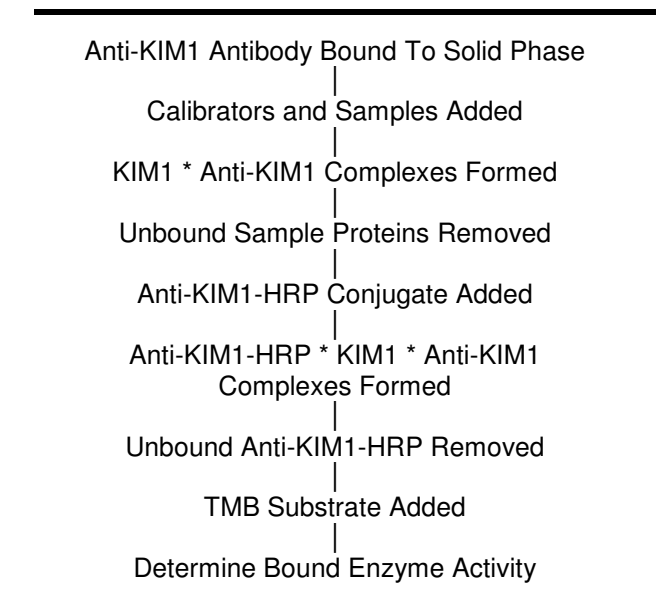
For Research Use Only.

PRODUCT INFORMATION**Rat Kidney Injury Molecule (KIM1) ELISA**
Cat. No. KT-622**INTENDED USE**

The Rat KIM1 ELISA is a highly sensitive two-site enzyme-linked immunoassay (ELISA) for the quantitative determination of KIM1 in Rat biological fluids. For research use only.

PRINCIPLE

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the KIM1 present in samples reacts with the anti-KIM1 antibodies which have been adsorbed to the surface of polystyrene microtiter wells. After the removal of unbound proteins by washing, anti-KIM1 antibodies conjugated with horseradish peroxidase (HRP) are added. These enzyme-labeled antibodies form complexes with the previously bound KIM1. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of KIM1 in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of KIM1 in the test sample. The quantity of KIM1 in the test sample can be interpolated from the calibration curve constructed from the calibrators, and corrected for sample dilution.

Figure 1.**COMPONENTS**

1. Diluent Concentrate
One bottle containing 50 mL of a 5X concentrated diluent running buffer.
2. Wash Solution Concentrate
One bottle containing 50 mL of a 20X concentrated wash solution.
3. Enzyme-Antibody Conjugate Concentrate
One vial containing 150 µL of a 100X concentrated affinity-purified anti-Rat KIM1 antibody conjugated with HRP in stabilizing buffer.

4. TMB Substrate Solution
One vial containing 12 mL of TMB and hydrogen peroxide in citric acid buffer at pH 3.3.
5. Stop Solution
One vial containing 12 mL of 0.3 M sulfuric acid. WARNING: Avoid contact with skin.
6. Microtiter Plate
Twelve removable eight-well strips in well holder frame. Wells are coated with affinity-purified anti-Rat KIM1.
7. Rat KIM1 Calibrator
One vial containing a lyophilized Rat KIM1 Calibrator.

MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes (2 μ L to 200 μ L) for making and dispensing dilutions
- Test tubes
- Microplate washer/aspirator
- Distilled or de-ionized H₂O
- Microplate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer
- Vortex mixer

PRECAUTIONS

1. Read the instructions carefully before beginning the assay.
2. This kit is for research use only.
3. Great care has been taken to ensure the quality and reliability of this product. However, it is possible that in certain cases, unusual results may be obtained due to high levels of interfering factors.
4. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.
5. Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.
6. Other precautions:
 - Do not interchange kit components from different lots.
 - Do not use kit components beyond the expiration date.
 - Protect reagents from direct sunlight.
 - Do not pipette by mouth.
 - Do not eat, drink, smoke or apply cosmetics where reagents are used.
 - Avoid all contact with the reagents by using gloves.
 - Stop solution contains diluted sulfuric acid. Irritation to eyes and skin is possible. Flush with water after contact.

REAGENT PREPARATION

1. Diluent Concentrate
The Diluent solution supplied is a 5X concentrate and must be diluted 1:5 with distilled or de-ionized water.
2. Wash Solution Concentrate
The Wash Solution supplied is a 20X concentrate and must be diluted 1:20 with distilled or de-ionized water. Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.
3. Enzyme-Antibody Conjugate Concentrate
Calculate the required amount of working conjugate solution for each microtiter plate test strip by adding 10 μ L Enzyme-Antibody Conjugate to 990 μ L of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.
4. TMB Substrate Solution
Ready to use as supplied.
5. Stop Solution
Ready to use as supplied.

6. Microtiter Plate

Ready to use as supplied. Unseal microtiter pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.

7. Rat KIM1 Calibrator

Add 0.5 mL of distilled or de-ionized water to the lyophilized Rat KIM1 Calibrator and mix gently until dissolved. The calibrator is now at a concentration of 500 ng/mL (the reconstituted calibrator should be aliquoted and frozen if future use is intended). Prepare the Rat KIM1 Calibrators immediately prior to use according to the table below. Mix well between each step. Avoid foaming.

Calibrator	Concentration (ng/mL)	Calibrator Volume added to 1X Diluent	Volume of 1X Diluent
7	10	200 μ L Rat KIM1 Calibrator	300 μ L
6	5	250 μ L Calibrator 7	250 μ L
5	2.5	250 μ L Calibrator 6	250 μ L
4	1.25	250 μ L Calibrator 5	250 μ L
3	0.625	250 μ L Calibrator 4	250 μ L
2	0.313	250 μ L Calibrator 3	250 μ L
1	0.156	250 μ L Calibrator 2	250 μ L
0	0		500 μ L

STORAGE AND STABILITY

1. Complete Kit

The expiration date for the kit is stated on the outer label. The recommended storage temperature is 4°C. **Note: See long term storage recommendations below for the Rat KIM1 Calibrator.**

2. Diluent

The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4°C.

3. Wash Solution

The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions can be stored at room temperature (RT, 16-25°C) or at 4°C.

4. Enzyme-Antibody Conjugate

Undiluted anti-KIM1-HRP conjugate should be stored at 4°C and diluted immediately prior to use. The working conjugate solution is stable for up to 8 hours.

5. TMB Substrate Solution

The TMB Substrate Solution should be stored at 4°C and is stable until the expiration date.

6. Stop Solution

The Stop Solution should be stored at 4°C and is stable until the expiration date.

7. Microtiter Plate

Anti-Rat KIM1 coated wells are stable until the expiration date, and should be stored at 4°C in the sealed foil pouch with desiccant pack.

8. Rat KIM1 Calibrator

The lyophilized Rat KIM1 Calibrator should be stored at 4°C or frozen until reconstituted. The reconstituted calibrator should be aliquoted and stored frozen. Avoid multiple freeze/thaw cycles. The working calibrator solutions should be prepared immediately prior to use and are stable for up to 1 hour.

INDICATIONS OF INSTABILITY

If the test is performing correctly, the results observed with the calibrator solutions should be within 20% of the expected values.

ASSAY PROTOCOL

Dilution of Samples

The assay for quantification of KIM1 in samples requires that each test sample be diluted before use. For a single step determination a dilution of 1/2 is appropriate for most serum/plasma samples. For absolute quantification of samples that yield results outside the range of the calibration curve, a lesser or greater dilution might be required. If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.

To prepare a 1/2 dilution of sample, transfer 150 μ L of sample to 150 μ L of 1X Diluent. This gives you a 1/2 dilution of your sample. Mix thoroughly.

Procedure

Bring all reagents to RT before use.

1. Pipette 100 μ L of
 - Calibrator 0 (0.0 ng/mL) in duplicate
 - Calibrator 1 (0.156 ng/mL) in duplicate
 - Calibrator 2 (0.313 ng/mL) in duplicate
 - Calibrator 3 (0.625 ng/mL) in duplicate
 - Calibrator 4 (1.25 ng/mL) in duplicate
 - Calibrator 5 (2.5 ng/mL) in duplicate
 - Calibrator 6 (5 ng/mL) in duplicate
 - Calibrator 7 (10 ng/mL) in duplicate
2. Pipette 100 μ L of diluted sample (in duplicate) into pre-designated wells.
3. Incubate the Microtiter Plate at 22°C (RT) for one hundred and twenty (120 \pm 2) minutes. Keep plate covered and level during incubation.
4. Following incubation, aspirate the contents of the wells.
5. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with diluted Wash Solution, invert the plate and pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual Wash Solution. Repeat three times for a total of four washes.
6. Pipette 100 μ L of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at 22°C (RT) for sixty (60 \pm 2) minutes. Keep plate covered in the dark and level during incubation.
7. Wash and blot the wells as described in Steps 4 and 5.
8. Pipette 100 μ L of TMB Substrate Solution into each well.
9. Incubate in the dark at RT for precisely ten (10) minutes.
10. After ten minutes, add 100 μ L of Stop Solution to each well.
11. Determine the absorbance at 450 nm of the contents of each well. Zero the plate reader to air.

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

RESULTS

1. Subtract the average background value from the test values for each sample.
2. Using the results observed for the calibrators construct a calibration curve. The appropriate curve fit is that of a four-parameter logistics curve, although a second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from calibration curve. Correct for sample dilution factor to arrive at KIM1 concentration in original sample.

PERFORMANCE CHARACTERISTICS

In accord with good laboratory practice, the assays for specific KIM1 require meticulous quality control. Each laboratory should use routine quality control procedures to establish inter- and intra-assay precision and performance characteristics.

LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or de-ionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.

FOR RESEARCH USE ONLY

KAMIYA BIOMEDICAL COMPANY

12779 Gateway Drive, Seattle, WA 98168
Tel: (206) 575-8068 Fax: (206) 575-8094
Email: LifeScience@k-assay.com
www.k-assay.com