

**KAMIYA BIOMEDICAL COMPANY**

# Mouse Amyloid Beta Peptide 1-40 (AB1-40) ELISA

**For the quantitative determination of mouse AB1-40 in  
serum, plasma and other biological fluids**

**Cat. No. KT-6382**

**For Research Use Only. Not for use in diagnostic procedures.**

**Product Information****Mouse Amyloid Beta Peptide 1-40 (AB1-40) ELISA****Cat. No. KT-6382****INTENDED USE**

The kit is a sandwich enzyme immunoassay for the *in vitro* quantitative measurement of mouse AB1-40 in serum, plasma and other biological fluids. For research use only. Not for use in diagnostic procedures.

**COMPONENTS**

<b>Reagents</b>	<b>Quantity</b>
Pre-coated, ready to use 96-well strip plate	1
Calibrator (lyophilized)	2
Calibrator Diluent	1 × 20 mL
Detection Reagent A	1 × 120 µL
Detection Reagent B	1 × 120 µL
Assay Diluent A (2X concentrate)	1 × 6 mL
Assay Diluent B (2X concentrate)	1 × 6 mL
TMB Substrate	1 × 9 mL
Stop Solution	1 × 6 mL
Wash Buffer (30X concentrate)	1 × 20 mL
Plate sealer for 96 wells	4

**MATERIALS REQUIRED BUT NOT SUPPLIED**

1. Microplate reader with  $450 \pm 10$  nm filter.
2. Precision single and multi-channel pipettes and disposable tips.
3. Eppendorf Tubes for diluting samples.
4. De-ionized or distilled water.
5. Absorbent paper for blotting the microtiter plate.
6. Container for Wash Solution.

**STORAGE**

All the reagents should be kept according to the labels on vials. The **Calibrator**, **Detection Reagent A**, **Detection Reagent B** and the **96-well strip plate** should be stored at -20°C upon being received. The unused strips should be kept in a sealed bag with the desiccant provided to minimize exposure to damp air. Opened test kits will remain stable until the expiration date shown, provided it is stored as prescribed above.

**PRINCIPLE**

The microtiter plate provided in this kit has been pre-coated with an antibody specific to AB1-40. A monoclonal antibody specific for mouse AB1-40 has been pre-coated onto a microplate. A competitive inhibition reaction is launched between biotin labeled mouse AB1-40 (calibrators and samples) with the pre-coated antibody specific for mouse AB1-40. After incubation the unbound conjugate is washed off. Next, avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. The amount of bound HRP conjugate is inversely proportional to the concentration of AB1-40

in the sample. After the addition of the substrate solution, the intensity of the color developed is inversely proportional to the concentration of AB1-40 in the sample.

## SAMPLE COLLECTION AND STORAGE

### Serum

Use a serum separator tube and allow samples to clot for two hours at room temperature or overnight at 4°C before centrifugation for 20 minutes at approximately 1,000 x g. Assay freshly prepared serum immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.

### Plasma

Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge samples for 15 minutes at 1,000 x g at 4°C within 30 minutes of collection. Remove plasma and assay immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.

### Other biological fluids

Centrifuge samples for 20 minutes at 1,000 x g. Remove particulates and assay immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid repeated freeze/thaw cycles.



### Note:


1. Samples to be used within 5 days may be stored at 4°C, otherwise samples must be stored at -20°C (≤1 month) or -80°C (≤2 months) to avoid loss of bioactivity and contamination.
2. When performing the assay slowly bring samples to room temperature.
3. Sample hemolysis will influence the result, so hemolytic specimen can not be detected.

## REAGENT PREPARATION

Bring all kit components and samples to room temperature (18-25°C) before use.

### Calibrator

Reconstitute the **Calibrator** with 0.5 mL of **Calibrator Diluent**, kept for 10 minutes at room temperature, shake gently (not to foam). The concentration of the calibrator in the stock solution is 1,000 pg/mL. Please prepare 5 tubes containing 0.6 mL Calibrator Diluent and produce a dilution series according to the picture shown below. Mix each tube thoroughly before the next transfer. Set up 5 points of diluted calibrator such as 1,000, 333.33, 111.11, 37.04, 12.35 pg/mL, and the last EP tubes with **Calibrator Diluent** is the blank as 0 pg/mL.

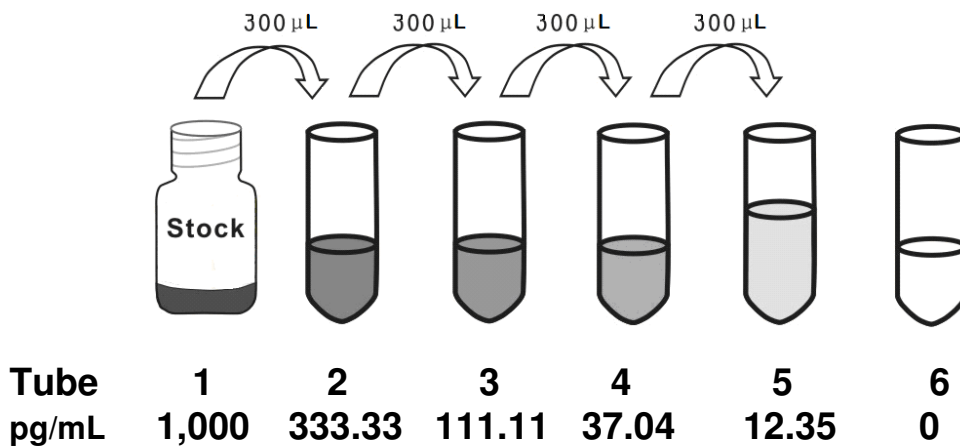

### Assay Diluent A and B

Dilute 6 mL of Assay Diluent A or B Concentrate (2X) with 6 mL of de-ionized or distilled water to prepare 12 mL of Assay Diluent A or B. **The prepared working dilution can't be frozen.**

### Detection Reagent A and B

Briefly spin or centrifuge the stock Detection Reagent A and Detection Reagent B before use. Dilute to the working concentration with working **Assay Diluent A or B**, respectively (1:100).

### Wash Solution

Dilute 20 mL of Wash Solution Concentrate (30X) with 580 mL of de-ionized or distilled water to prepare 600 mL of Wash Solution (1X).

### TMB substrate

Aspirate the needed dosage of the solution with sterilized tips and do not dump the residual solution into the vial again.



### Note:

1. Making serial dilution in the wells directly is not permitted.
2. Prepare calibrator within 15 minutes before assay. Please do not dissolve the reagents at 37°C directly.
3. Please carefully reconstitute Calibrators or working Detection Reagent A and B according to the instruction, and avoid foaming and mix gently until the crystals have completely dissolved. To minimize imprecision caused by pipetting, use small volumes and ensure that pipettors are calibrated. It is recommended to draw more than 10 µL for pipetting.
4. The reconstituted Calibrators, Detection Reagent A and Detection Reagent B can be **used only once**.
5. If crystals have formed in the Wash Solution concentrate (30X), warm to room temperature and mix gently until the crystals have completely dissolved.
6. Distilled water is recommended to be used to make the dilution for reagents or samples. Contaminated water or container for reagent preparation will influence the detection result.

## ASSAY PROCEDURE

1. Determine wells for diluted calibrator, blank and sample. Prepare 5 wells for calibrator, 1 well for blank. Add 50 µL each of dilutions of calibrator (read Reagent Preparation), blank and samples into the appropriate wells. Then add 50 µL of Detection Reagent A to each well immediately. Shake plate gently. Cover with the Plate sealer. Incubate for 1 hour at 37°C. Detection Reagent A may appear cloudy. Warm to room temperature and mix gently until solution appears uniform.
2. Aspirate the solution and wash with 350 µL of 1X Wash Solution to each well using a squirt bottle, multi-channel pipette, manifold dispenser or autowasher, and let it sit for 1~2 minutes. Remove the remaining liquid from all wells completely by snapping the plate onto absorbent paper. Repeat 3 times. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against absorbent paper.
3. Add 100 µL of **Detection Reagent B** working solution to each well. Incubate for 30 minutes at 37°C after covering it with the Plate sealer.
4. Repeat the aspiration/wash process for five times as conducted in step 2.
5. Add 90 µL of **Substrate Solution** to each well. Cover with a new Plate sealer. Incubate for 15 - 25 minutes at 37°C (Don't exceed 30 minutes). Protect from light. The liquid will turn blue by the addition of Substrate Solution.
6. Add 50 µL of **Stop Solution** to each well. The liquid will turn yellow by the addition of Stop solution. Mix the liquid by tapping the side of the plate. If color change does not appear uniform, gently tap the plate to ensure thorough mixing.
7. Remove any drop of water and fingerprint on the bottom of the plate and confirm there is no bubble on the surface of the liquid. Then, run the microplate reader and conduct measurement at 450 nm immediately.

**Note:**

1. **Assay preparation:** Keep appropriate numbers of strips for 1 experiment and remove extra strips from microtiter plate. Removed strips should be resealed and stored at -20°C until the kit expiration date.
2. **Samples or reagents addition: Please use the freshly prepared Calibrator.** Please carefully add samples to wells and mix gently to avoid foaming. Do not touch the well wall if possible. For each step in the procedure, total dispensing time for addition of reagents or samples to the assay plate should not exceed 10 minutes. This will ensure equal elapsed time for each pipetting step, without interruption. Duplication of all calibrators and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each calibrator level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
3. **Incubation:** To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
4. **Washing:** The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading.
5. **Controlling of reaction time:** Observe the change of color after adding **TMB Substrate** (e.g. observation once every 10 minutes), if the color is too deep, add **Stop Solution** in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
6. **TMB Substrate** is easily contaminated. Please protect it from light.

## CALCULATION OF RESULTS

This assay employs a competitive inhibition enzyme immunoassay technique, so there is an inverse correlation between AB1-40 concentration in the sample and the assay signal intensity. Low levels of AB1-40 result in a high OD value, while high concentrations result in a low signal.

Average the duplicate readings for each calibrator, control, and samples. Create a calibration curve on log-log graph paper, with AB1-40 concentration on the y-axis and absorbance on the x-axis. Draw the best fit straight line through the calibrator points and it can be determined by regression analysis. Using some plot software is also recommended. If samples have been diluted, the concentration read from the calibration curve must be multiplied by the dilution factor.

## PERFORMANCE

### Detection Range

12.35 - 1,000 pg/mL.

The calibration curve concentrations used for the ELISA's were 1,000, 333.33, 111.11, 37.04, 12.35 pg/mL.

### Sensitivity

The minimum detectable dose of mouse AB1-40 is typically less than 5.73 pg/mL. The sensitivity of this assay, or Lower Limit of Detection (LLD) was defined as the lowest protein concentration that could be differentiated from zero. It was determined the mean O.D. Value of 20 replicates of the zero calibrator plus three standard deviations.

### Specificity

This assay has high sensitivity and excellent specificity for detection of mouse AB1-40. No significant cross-reactivity or interference between mouse AB1-40 and analogues was observed.

Note: Limited by current skills and knowledge, it is impossible for us to complete the cross-reactivity detection between mouse AB1-40 and all the other analogues, therefore, cross reaction may still exist.

**IMPORTANT NOTES**

1. The final experimental results will be closely related to operation skills of the end users and the experimental environments. Please make sure that sufficient samples are available.
2. Kits from different batches may be a little different in detection range, sensitivity and color developing time. Please perform the experiment exactly according to the instruction attached in kit while electronic ones from our website ([www.k-assay.com](http://www.k-assay.com)) is only for information.
3. Do not mix or substitute reagents from one kit lot to another. Use only the reagents supplied by manufacturer.
4. Protect all reagents from strong light during storage and incubation. All the bottle caps of reagents should be covered tightly to prevent the evaporation and contamination of microorganism.
5. There may be some foggy substance in the wells when the plate is opened at the first time. It will not have any effect on the final assay results. Do not remove microtiter plate from the storage bag until needed.
6. Wrong operations during the reagents preparation and loading, as well as incorrect parameter setting for the plate reader may lead to incorrect results. A microplate plate reader with a bandwidth of 10 nm or less and an optical density range of 0-3 O.D. or greater at  $450 \pm 10$  nm wavelength is acceptable for use in absorbance measurement. Please read the instruction carefully and adjust the instrument prior to the experiment.
7. Even the same operator might get different results in two separate experiments. In order to get better reproducible results, the operation of every step in the assay should be controlled. Furthermore, a preliminary experiment before assay for each batch is recommended.
8. Each kit has been strictly passed Q.C. test. However, results from end users might be inconsistent with our in-house data due to some unexpected transportation conditions or different lab equipments. Intra-assay variance among kits from different batches might arise from above factors, too.
9. Kits from different manufacturers for the same item might produce different results, since we haven't compared our products with other manufacturers.
10. The Stop Solution suggested for use with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material.

**FOR RESEARCH USE ONLY. Not for use in diagnostic procedures.**

**KAMIYA BIOMEDICAL COMPANY**

12779 Gateway Drive, Seattle, WA 98168  
Tel: (206) 575-8068 Fax: (206) 575-8094  
Email: [LifeScience@k-assay.com](mailto:LifeScience@k-assay.com)  
[www.k-assay.com](http://www.k-assay.com)