

KAMIYA BIOMEDICAL COMPANY

Guinea Pig IgM ELISA

For the quantitative determination of IgM in guinea pig serum or plasma

Cat. No. KT-656

For Research Use Only.

PRODUCT INFORMATION

Guinea Pig IgM ELISA **Cat. No. KT-656**

PRODUCT

The **K-ASSAY®** Guinea Pig IgM ELISA is an enzyme immunoassay for the quantitative determination of IgM in guinea pig serum or plasma. For research use only.

INTRODUCTION

The guinea pig IgM ELISA kit is designed for measurement of IgM in guinea pig serum or plasma. The assay uses goat anti-guinea pig IgM for solid phase (microtiter wells) immobilization and horseradish peroxidase (HRP) conjugated goat anti-guinea pig IgM antibodies for detection.

PRINCIPLE

Test samples are diluted and incubated in the microtiter wells for 45 minutes alongside prepared guinea pig IgM calibrators. The microtiter wells are subsequently washed and HRP conjugate is added and incubated for 45 minutes. IgM molecules are thus sandwiched between the immobilization and detection antibodies. The wells are then washed to remove unbound HRP-labeled antibodies and TMB Reagent is added and incubated for 20 minutes at room temperature. This results in the development of a blue color. Color development is stopped by the addition of Stop Solution, changing the color to yellow, and optical density is measured spectrophotometrically at 450 nm. The concentration of IgM is proportional to the optical density of the test sample and is derived from a calibration curve.

COMPONENTS

- Anti guinea pig IgM coated 96-well plate (12 strips of 8 wells)
- Reference calibrator (lyophilized)
- 10X Diluent, 25 mL
- HRP Conjugate Reagent, 11 mL
- 20X Wash Solution, 50 mL
- TMB Reagent (One-step), 11 mL
- Stop Solution (1N HCl), 11 mL

MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes and tips
- Distilled or de-ionized water
- Vortex mixer
- Absorbent paper or paper towels
- Graph paper (PC graphing software is optional)
- Polypropylene or glass tubes
- Plate reader with an optical density range of 0-4 at 450 nm.
- Micro-Plate incubator/shaker mixing speed of ~150 rpm
- Plate washer

DILUENT PREPARATION

The diluent is provided as a 10X stock. Prior to use estimate the final volume of diluent required for your assay and dilute one (1) volume of the 10x stock with nine (9) volumes of distilled or de-ionized water.

WASH SOLUTION PREPARATION

The wash solution is provided as a 20X stock. Prior to use dilute the contents of the bottle (50 mL) with 950 mL of distilled or de-ionized water.

CALIBRATOR PREPARATION

1. The IgM calibrator is provided as a lyophilized stock. Reconstitute as detailed on the vial label (***the reconstituted calibrator is stable at 4°C for one week but should be aliquoted and frozen at -20°C after reconstitution if future use is intended***).
2. Label 7 polypropylene or glass tubes as 250, 125, 62.5, 31.25, 15.63, 7.81, and 3.91 ng/mL.
3. Into the tube labeled 250 ng/mL, pipette the volume of diluent detailed on the IgM calibrator vial label. Then add the indicated volume of IgM calibrator and mix gently. This provides the 250 ng/mL calibrator.
4. Dispense 250 μ L of diluent into the tubes labeled 125, 62.5, 31.25, 15.63, 7.81, and 3.91 ng/mL.
5. Prepare a 125 ng/mL calibrator by diluting and mixing 250 μ L of the 250 ng/mL calibrator with 250 μ L of diluent in the tube labeled 125 ng/mL.
6. Similarly prepare the 62.5, 31.25, 15.63, 7.81, and 3.91 ng/mL calibrators by serial dilution.

SAMPLE PREPARATION

General Note: IgM is typically present in guinea pig serum at concentrations of ~0.3 mg/mL. In order to obtain values within range of the calibration curve, we suggest that samples initially be diluted 5,000 fold using the following procedure for each sample to be tested:

1. Dispense 198 μ L and 294 μ L of 1X diluent into separate tubes.
2. Pipette and mix 2 μ L of the serum/plasma sample into the tube containing 198 μ L of diluent. This provides a 100 fold diluted sample.
3. Mix 6 μ L of the 100 fold diluted sample with the 294 μ L of diluent in the second tube. This provides a 5,000 fold dilution of the sample.
4. Repeat this procedure for each sample to be tested.

Tissue extracts and tissue culture medium will likely have lower IgM levels than those found in serum. Optimal dilutions of such samples should be determined empirically.

ASSAY PROCEDURE

1. Secure the desired number of coated wells in the holder.
2. Dispense 100 μ L of calibrators and diluted samples into the wells (we recommend that samples be tested in duplicate).
3. Incubate on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 45 minutes.
4. Aspirate the contents of the microtiter wells and wash the wells 5 times with 1x wash solution using a plate washer (400 μ L/well). The entire wash procedure should be performed as quickly as possible.
5. Strike the wells sharply onto absorbent paper or paper towels to remove all residual wash buffer.
6. Add 100 μ L of enzyme conjugate reagent into each well.
7. Incubate on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 45 minutes.
8. Wash as detailed in 4 to 5 above.
9. Dispense 100 μ L of TMB Reagent into each well.
10. Gently mix on an orbital micro-plate shaker at 100-150 rpm at room temperature (18-25°C) for 20 minutes.
11. Stop the reaction by adding 100 μ L of Stop Solution to each well.
12. Gently mix. *It is important to make sure that all the blue color changes to yellow.*
13. Read the optical density at 450 nm with a microtiter plate reader within 5 minutes.


CALCULATION OF RESULTS

1. Calculate the average absorbance values (A_{450}) for each set of reference calibrators and samples.
2. Construct a calibration curve by plotting the mean absorbance obtained from each reference calibrator against its concentration in ng/mL on linear graph paper, with absorbance values on the vertical or Y-axis and concentrations on the horizontal or X-axis.
3. Using the mean absorbance value for each sample, determine the corresponding concentration of IgM in ng/mL from the calibration curve.
4. Multiply the derived concentrations by the dilution factor to determine the actual concentration of IgM in the sample.
5. PC graphing software may be used for the above steps.
6. If the OD_{450} values of the sample fall outside the calibration curve, samples should be diluted appropriately and re-tested.

TYPICAL CALIBRATION CURVE

A typical calibration curve with optical density readings at 450 nm on the Y axis against IgM concentrations on the X axis is shown below. This curve is for the purpose of illustration only and should not be used to calculate unknowns. Each user should obtain his or her data and calibration curve in each experiment.

IgM (ng/mL)	Absorbance (450 nm)
250	2.5040
125	1.5250
62.5	0.9120
31.3	0.4920
15.63	0.3040
7.81	0.1900
3.91	0.1210

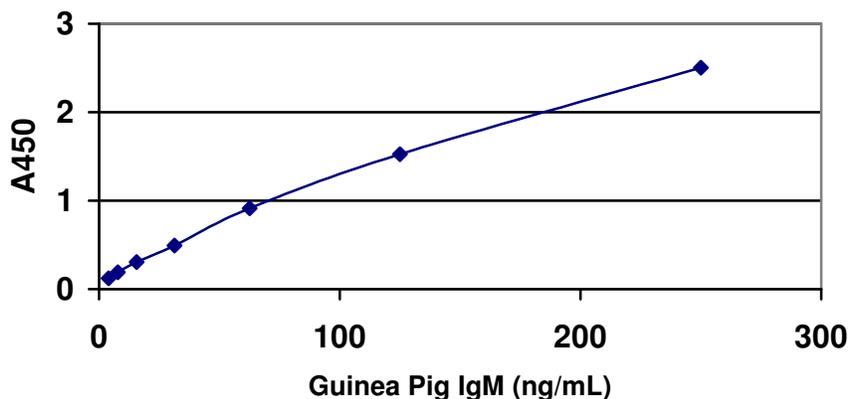

STORAGE

The test kit will remain stable until the expiration date provided that the components are stored at 4°C. The microtiter plate should be kept in a sealed bag with desiccant to minimize exposure to damp air.

GENERAL INSTRUCTIONS

1. Please read and understand the instructions thoroughly before using the kit.
2. All reagents should be allowed to reach room temperature (18-25°C) before use.
3. Optimum results are achieved if, at each step, reagents are pipetted into the wells of the microtiter plate within 5 minutes.

LIMITATIONS OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of and in accordance with the instructions detailed above.
2. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.

FOR RESEARCH USE ONLY

KAMIYA BIOMEDICAL COMPANY

12779 Gateway Drive, Seattle, WA 98168
 Tel: (206) 575-8068 Fax: (206) 575-8094
 Email: LifeScience@k-assay.com
 www.k-assay.com